



User Manual

Mandala Wand



Please read the manual before using the tools.

Version 5.2
Updated July 2, 2026

1. General Notes

1.1 Product Feature and Intend Use

Mandala Wand is a sleek, portable wellness device designed to restore balance and enhance well-being using advanced light therapy and scalar wave technology. With its powerful quantum and bifilar coils, it helps harmonize your body's energy flow, promoting relaxation and vitality.

Pair it with the Mandala Wand app, which allows you to upload and create custom presets for personalized energy balancing. The app also includes built-in preset programs for immediate use. Additionally, by integrating Mandala Scan, you can analyze your body's energy state and identify the best coherence preset to upload for optimal results.

Your Mandala Wand has been quality-tested and is ready to support your wellness journey anytime, anywhere.

1.2 Product details

- Light Frequency Modulator for Therapy
- Uses electronic control circuit
- Uses RGB LED, UV-C LED, IR-LASER, VIOLET – LASER, PCB Coil and PCB ANT as output
- Support Bluetooth for mobile control
- Uses 3.5 mm Jack Output to listen to music through the device

1.3 Safety Warnings

- Check all accessories upon unboxing.
- Read and follow all safety instructions in the manual.
- Do not turn the device on by looking directly at the device because it contains lasers.
- Do not turn on near water or in wet areas.
- Do not take apart if you are not a professional technician.
- Avoid placing near heat sources such as stoves.
- Use only the power supply specified in the manual.

1.4 Warranty

- Your Mandala Wand includes a one-year warranty from the date of purchase.
- If you experience any issues, contact customer support immediately.
- The warranty does not cover damage due to:
 - Misuse or neglect
 - Unauthorized modifications or disassembly
 - Customs/import fees or indirect damage
- If repairs are needed due to misuse, all costs (including shipping) will be the user's responsibility.

1.5 Important Precautions

- Do not use the device if it appears damaged. Contact support.
- Ensure the correct voltage for your country before use.
- Charge for at least 1 and half hours before the first use.
- Do not disassemble or modify the device.
- Keep dry – avoid rain, high humidity, heat, and dust.

1.6 Maintenance

- If unused for a while, charge the device for 2 hours before use.
- Use soft cloth or water-based glass cleaner for maintenance.
- Keep the device clean to ensure proper performance.
- Keep in dry areas, avoid direct sunlight
- If you have any problems, please contact customer service

1.7 Getting started

- To turn on the device, press the button once. The device will start working according to the basic functions.
- Connect to your mobile phone via Bluetooth: MW2-XXXX
- Select the desired treatment program via the application.

- To turn off the device, press the button twice in a quick double-press. The device will stop working immediately.

1.8 Advanced Setting

- There is a Mode button with 4 modes to select: Frequency, Basic Program/Preset, Bluetooth, and Audio.
- There is a Select button to select Basic program and preset that we have uploaded to the device in all 5 channels.
- There is a Start button
 - Start/stop the operation of Basic program and preset.
- In the application:
 - There are 16 groups of Basic Programs to choose from.
 - The Frequency mode function has 41 Library groups to choose from.
 - The Frequency mode has a function to enter the frequency number users want from 1-25MHz and can adjust the LED light.
 - There is a function that can be set via Mandala scan.
 - There is a Preset function that can create Presets that we want as our own and create and upload to store in Mandala Wand up to 5 Presets.
 - To use the application, Bluetooth mode needs to always be on.
 - There is a function that can set the language in a total of 11 languages.
- After opening the function, the machine will automatically shut itself off after 45 minutes of operation per usage.

1.9 Charging System

- When disconnecting the charger - normal use: there will be 4 LED light colors indicating the battery status: green, green, orange, red, representing the battery percentage 100% 75% 50% 25% in order

- When the battery is almost empty, only the red LED will blink to notify the battery is less than 10%
- When charging the battery, the device will not work. The LED will blink according to the battery percentage

1.10 Product Includes

- Mandala Wand
- USB Type-C cable
- Mandala Case
- 1x Symbol Sphere
- Crystal holder made of brass




2. Technical Data

2.1 Specifications

Product name:	Mandala Wand
Case Color:	Gold
Material:	ABS
Input Voltage:	5V DC
Charging Interface:	USB Type-C
Battery Type:	Lithium polymer 3.7 V/ 2000mA
Operation:	Bluetooth 5.0
Product Size:	57 x 175 x 32 mm (W x H x D)
Gross Weight:	142g
Main Frequency:	1 Hz – 25MHz

2.2 Overview and Function

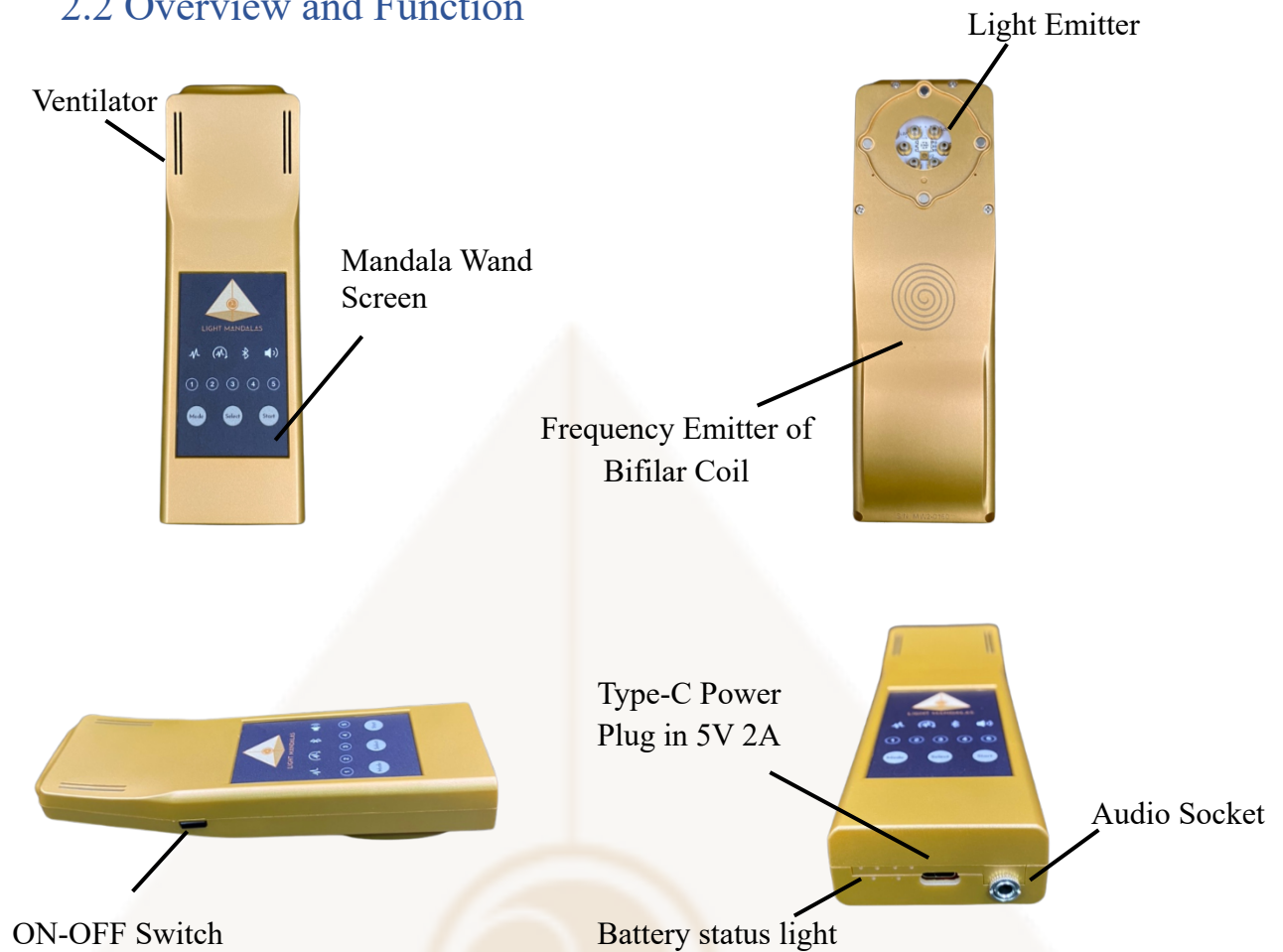

2.3 Crystal Sphere Holder Setup

The holder is designed to fit perfectly on the Mandala Wand. Every holder has a mark on one of the four tips. The mark serves as a guide for the user to align and place the Mandala Wand's top magnet.



2.4 Charging Setup

1. To charge the Mandala Wand, connect the USB-C cable to the device and a power adapter, then plug it into a power source.
2. Charging takes about 1 and half hours from an empty battery and lasts over 4 hours. (*Battery capacity can be reduced.*)
3. When the battery is low, the status light will be entirely red. When the battery is fully charged, the status light will display red, orange, and green.



To protect the battery life span, please don't overcharge it.

2.5 Mandala Wand Screen Function

When the device is switched on, wait a few seconds for the initialization process to be completed.

Preset List Mode:

Users can select created custom presets or from a Basic Program which has been uploaded to the Mandala Wand. This allows playback of presets without needing a phone or tablet in the future.

Frequency List Mode:


This mode allows users to select and play specific frequencies. Users can also upload frequencies from Frequency Library to the Mandala Wand, enabling playback without a phone or tablet in the future.

Bluetooth Mode:

Before pairing the Bluetooth with the app, please switch the mode to "Bluetooth Mode." This mode needs to be selected every time the application is being used, otherwise the device will disconnect.

Audio Mode:

Switch to Audio Mode before playing healing music. Users can control the Mandala Wand via the app and play music from the tablet. The lasers and scalar waves will pulse in sync with the music.

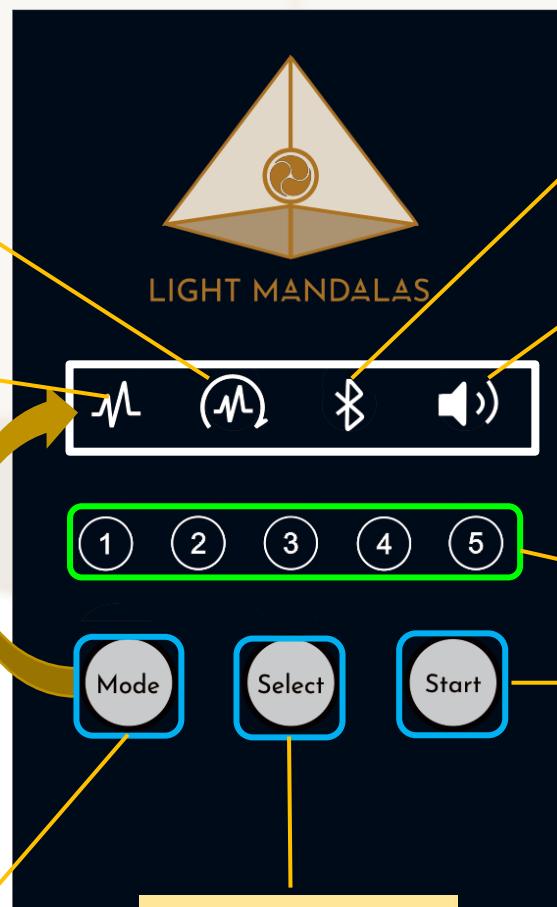

5 tools of frequency

Start is to start playing the presets. (It is used for Preset List mode)

Mode is for selecting 4 modes above.


Select is for selecting frequency presets from slot 1 to slot 5.

3. Mandala Wand Application


3.1 Application Setup

The Mandala Wand app is available on iOS and Android. Download it from the App Store for iOS or Google Play Store for Android. For users in China, please download it from <https://light-mandalas.com/application/>

App Store (iOS)

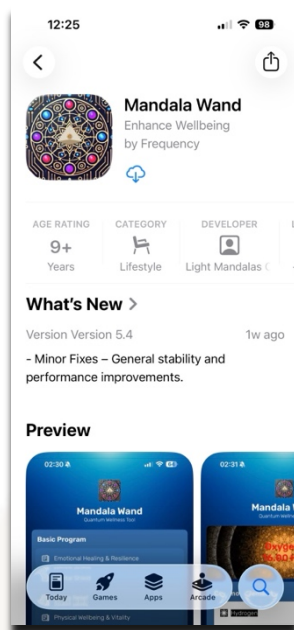


Play Store (Android)

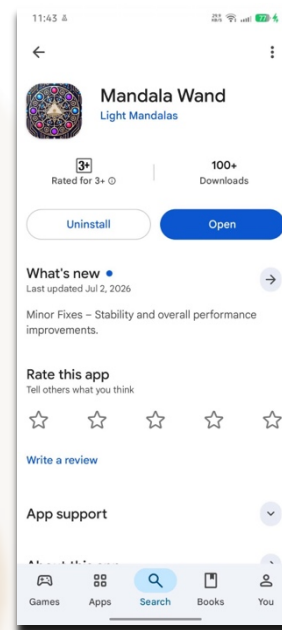


3.2 Pairing Bluetooth with Application

Select "Bluetooth" Mode  on your Mandala Wand. Tap the "Settings" icon on the bottom right, then select "Setup Mandala Wand." The application will search for nearby devices. Choose the displayed serial number to connect via Bluetooth.

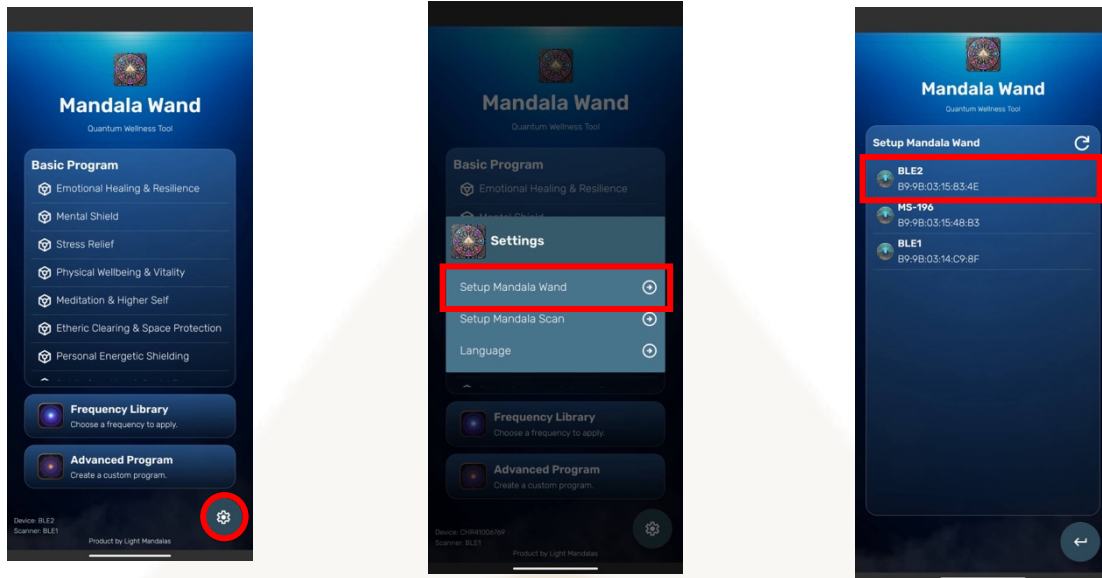

3.3 Basic Program

After setting up, select your preferred program under "Basic Program" as shown in the picture below. You can either play the program or upload it to a slot (1-5) for future use.

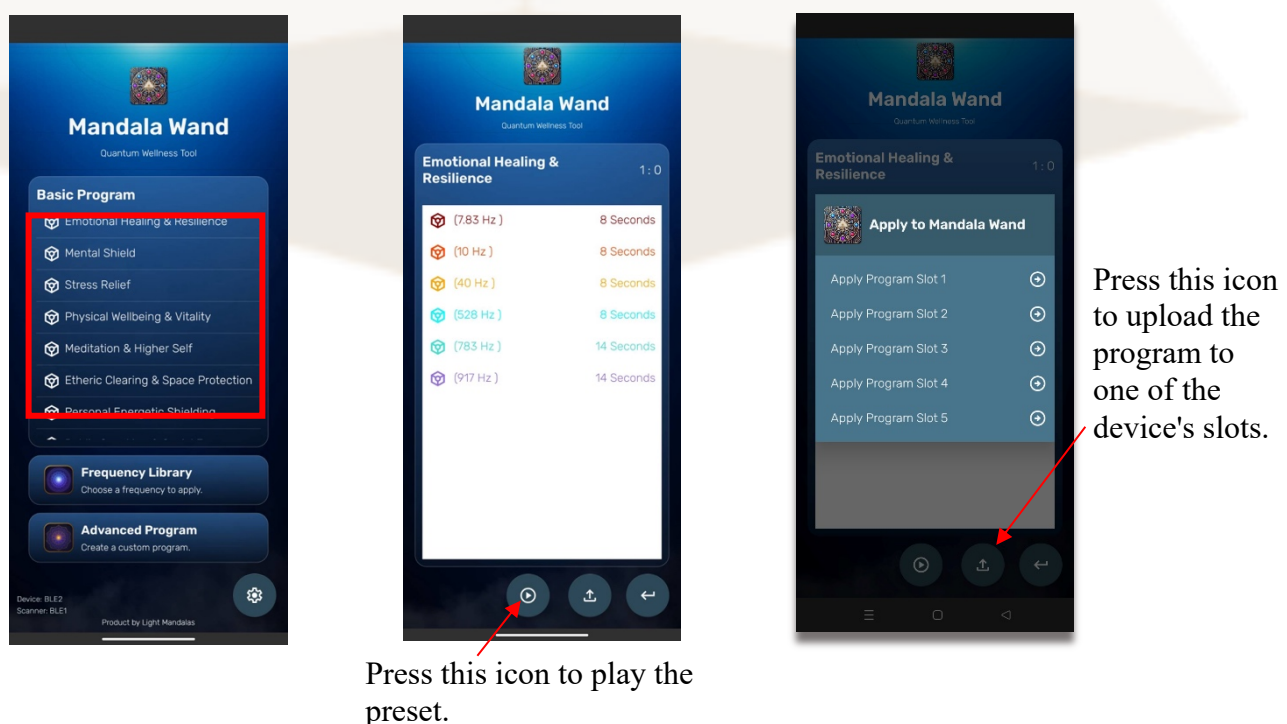


3.4 Selecting Mode and Slots

Programs uploaded to the slots can be played without using the application. To select the different programs uploaded to the slots, please do the following:

2. Select "Preset" mode.

1. Press "Mode" to select the mode you want to pick. To select from the Basic Program which has been uploaded to a slot, select "Preset" mode.

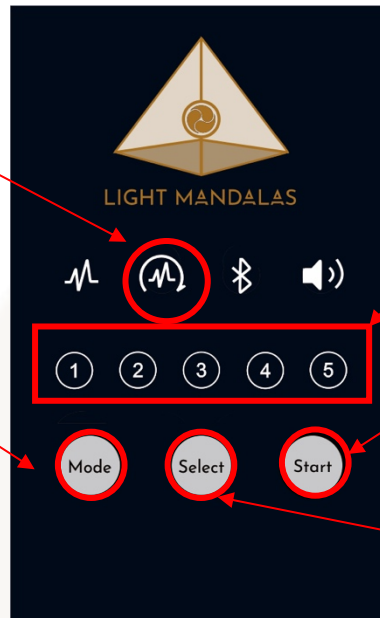


4. Continue to press "Select". The light will switch between the numbers.



5. Press "Start" to start playing the program.

3. Press "Select" to pick the slot number.

3.4.1 Frequency Library

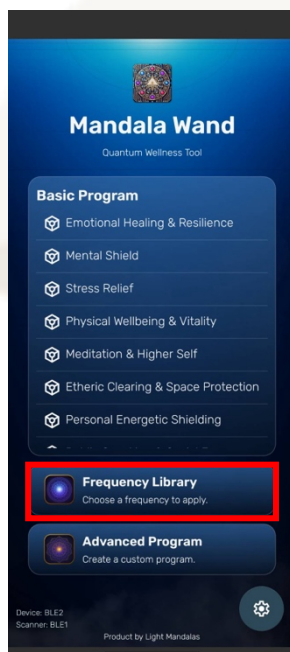

You can pick from Frequency Library.

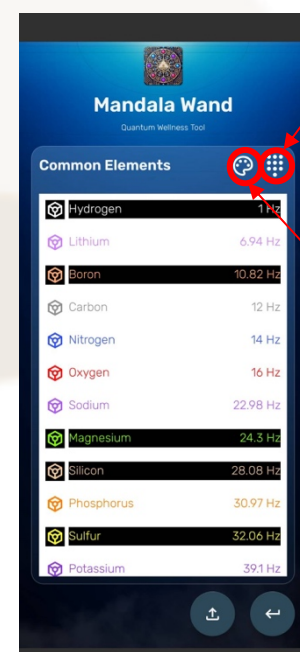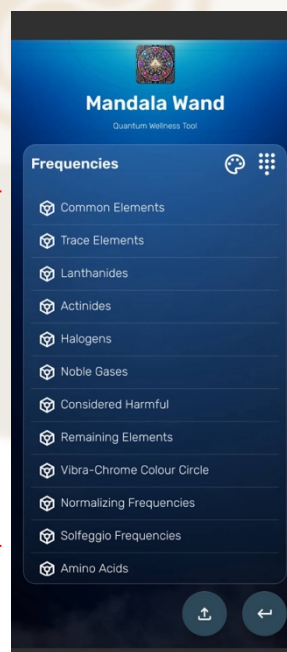

Insert your own frequency.

Pick a color for your custom frequency.

3.4.2 Insert Frequency and Frequency Color



Insert Frequency



Press this icon to insert your own frequency and "Confirm."

Press the color palette icon to pick your desired color. Use the slider to adjust the shade and then press "Confirm."




Frequency Color

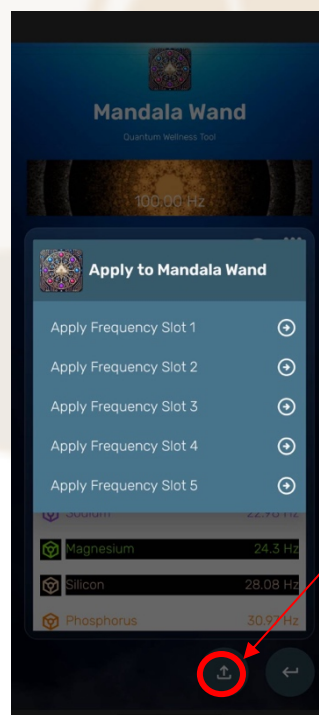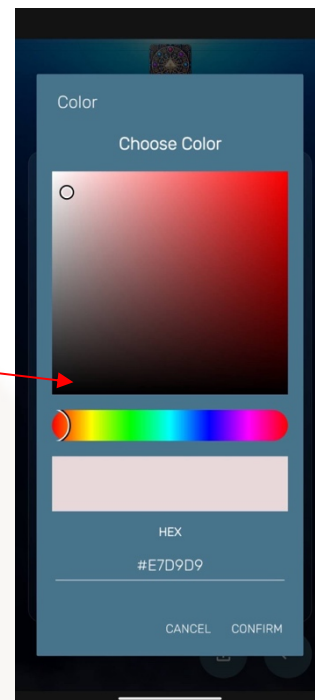



Press this icon to upload your preferred Frequency to the Frequency slot. After upload, you can play the Frequency offline.

3.5 Advanced Program


Advanced Program allows users to use the Mandala Wand with the Mandala Scan, as well as create their own custom programs and pick colors.

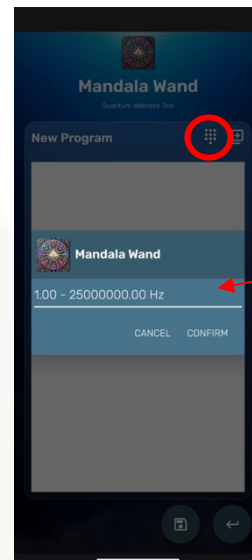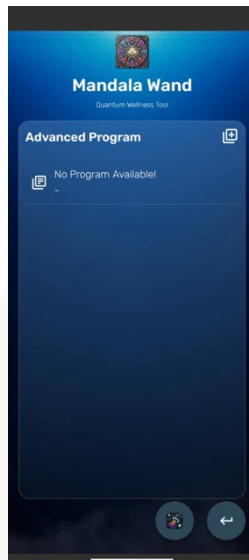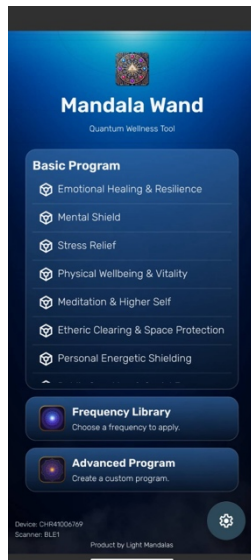

Insert your own frequency to put in your custom program.

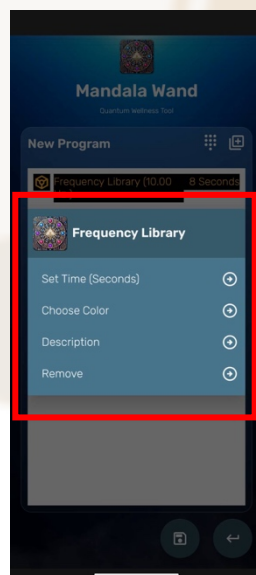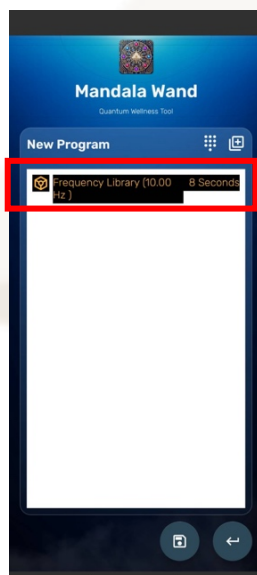

By pressing onto their custom frequency, users can set time, choose colors, write a description or remove their personal frequency.

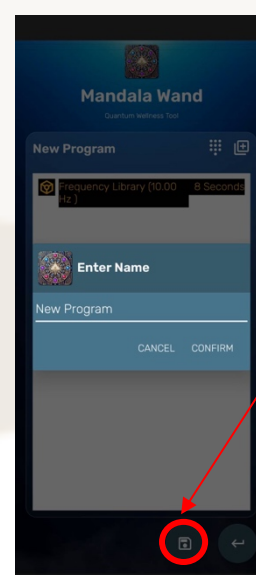

Press this icon to save the program created. Name the program and then press "Confirm."







3.5.1 Mandala Scan

To use Mandala Scan, first pair it via Bluetooth. Tap the "Settings" icon, select "Setup Mandala Scan," and choose the displayed serial number to connect.



3.5.2 Smart Scan

The Smart Scan function scans you for the best frequency for your body. The results are frequencies which your body needs in the current moment. You can play and upload the results to your Mandala Wand.

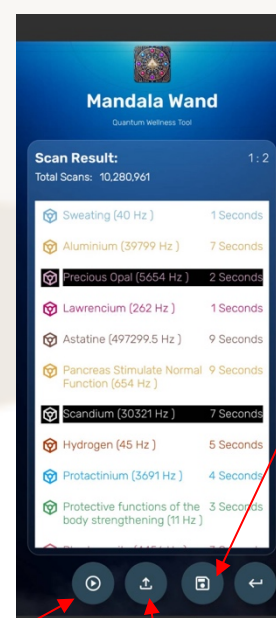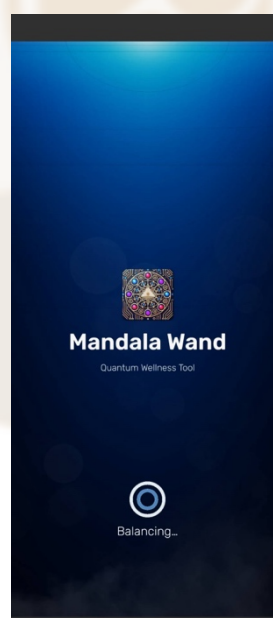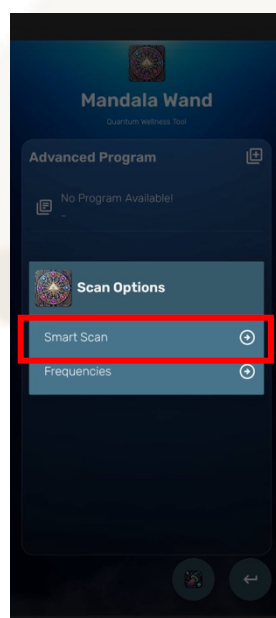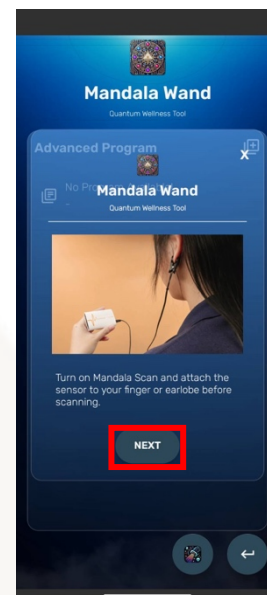
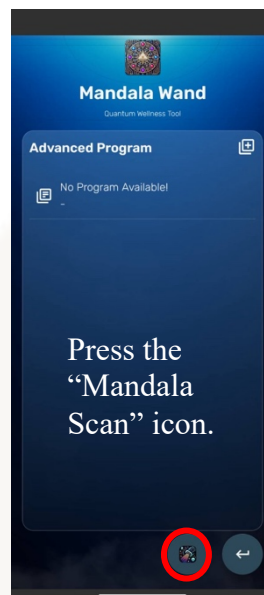
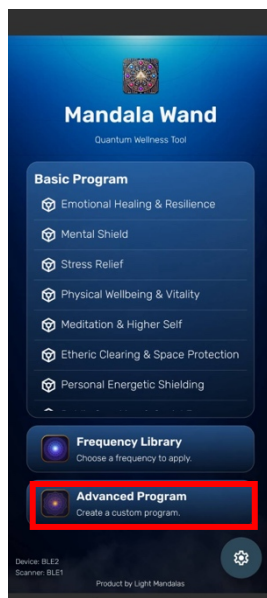


Press this icon to save the result into a custom program. Here, you can name your program, and it will appear in the library page.


Press this icon to play the scan results.

Press this icon to upload the scanning result to the Mandala Wand slots.


3.5.3 Frequencies Scan

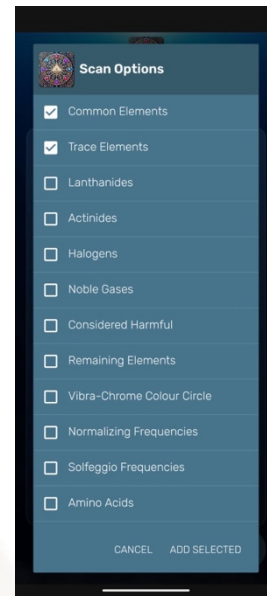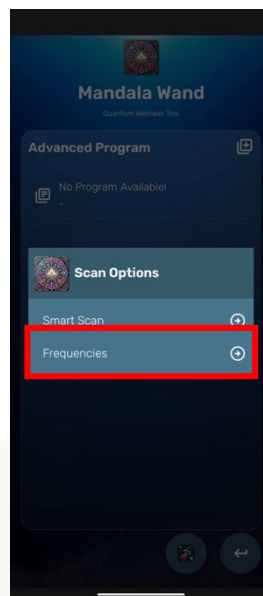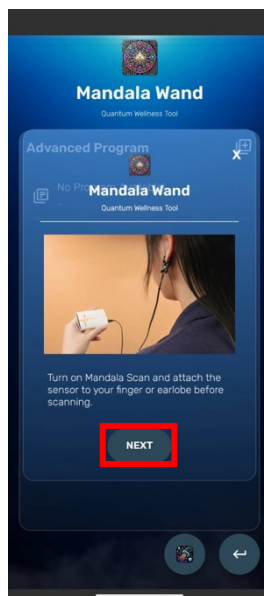

Select the frequency that you want to scan. You can add multiple categories. Press “Add Selected.”

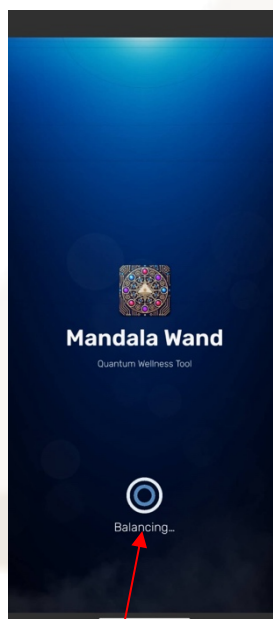

During scanning, the color will change from Red, to Blue and lastly to Green.


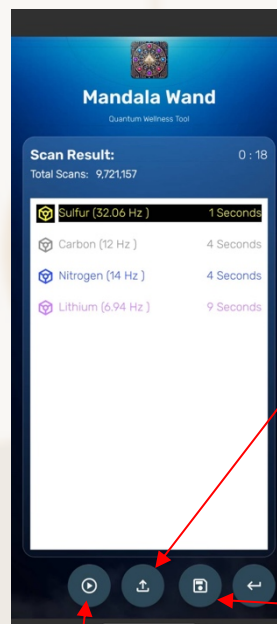


Press this icon to play the scanning result.

Press this icon to upload the scanning result to the Mandala Wand slots.

Press this icon to save the result into a custom program. Here, you can name your program, and it will appear in the library page.

3.5.4 Custom Program

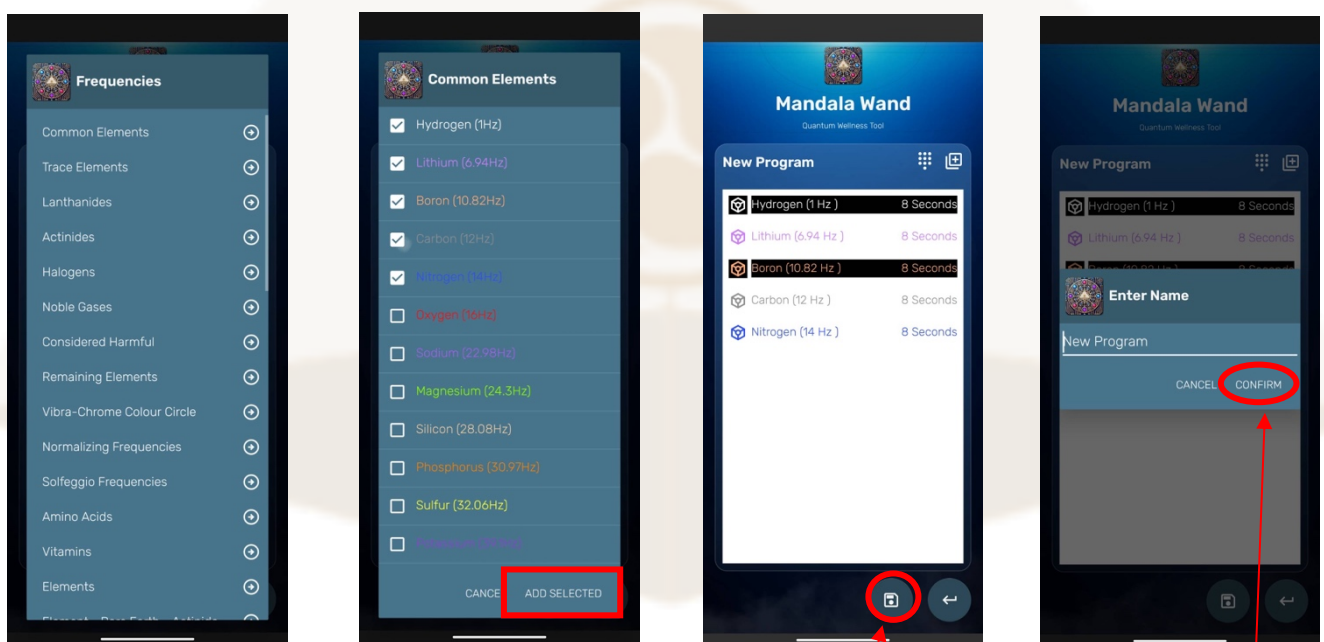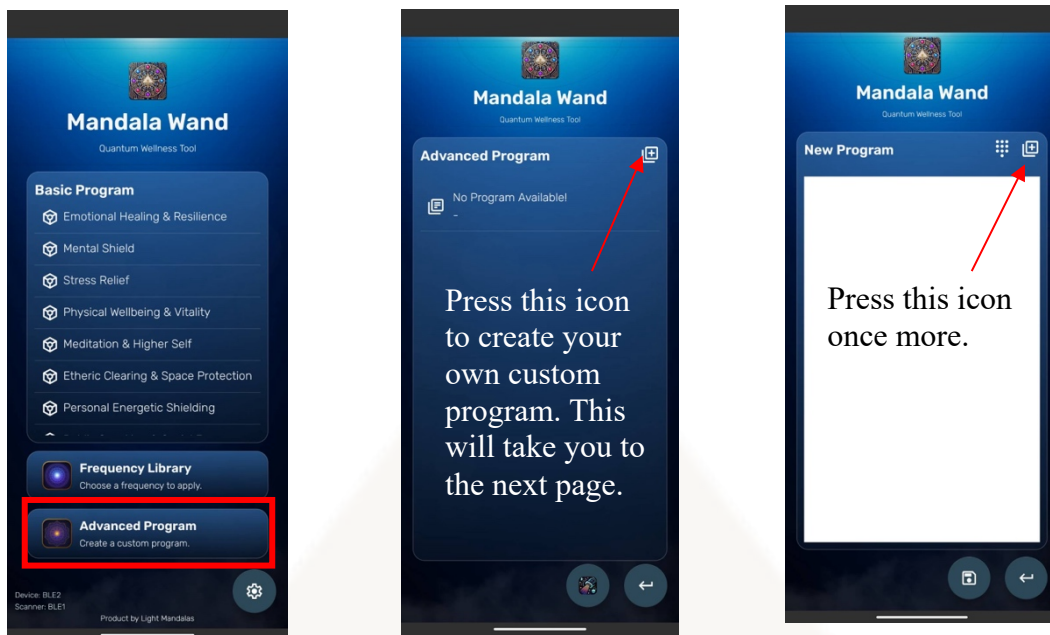


This is the list of Frequencies. Users can select their preferred Frequency and press “Add Selected.”

Press this icon to save your custom program.

Enter a name for your custom program and then press “Confirm.”

3.6 Audio Frequency

Switch to **Audio Mode** on the device's screen. Connect your headphones or speaker to the Mandala Wand, then pair the Mandala Wand Audio Bluetooth in your Bluetooth settings. Press “XFW-BT” to connect. Once you hear the notification sound through your headphones, you can start playing healing music from any music platform. Be sure to turn the volume on high to be able to see the pulsing colors.

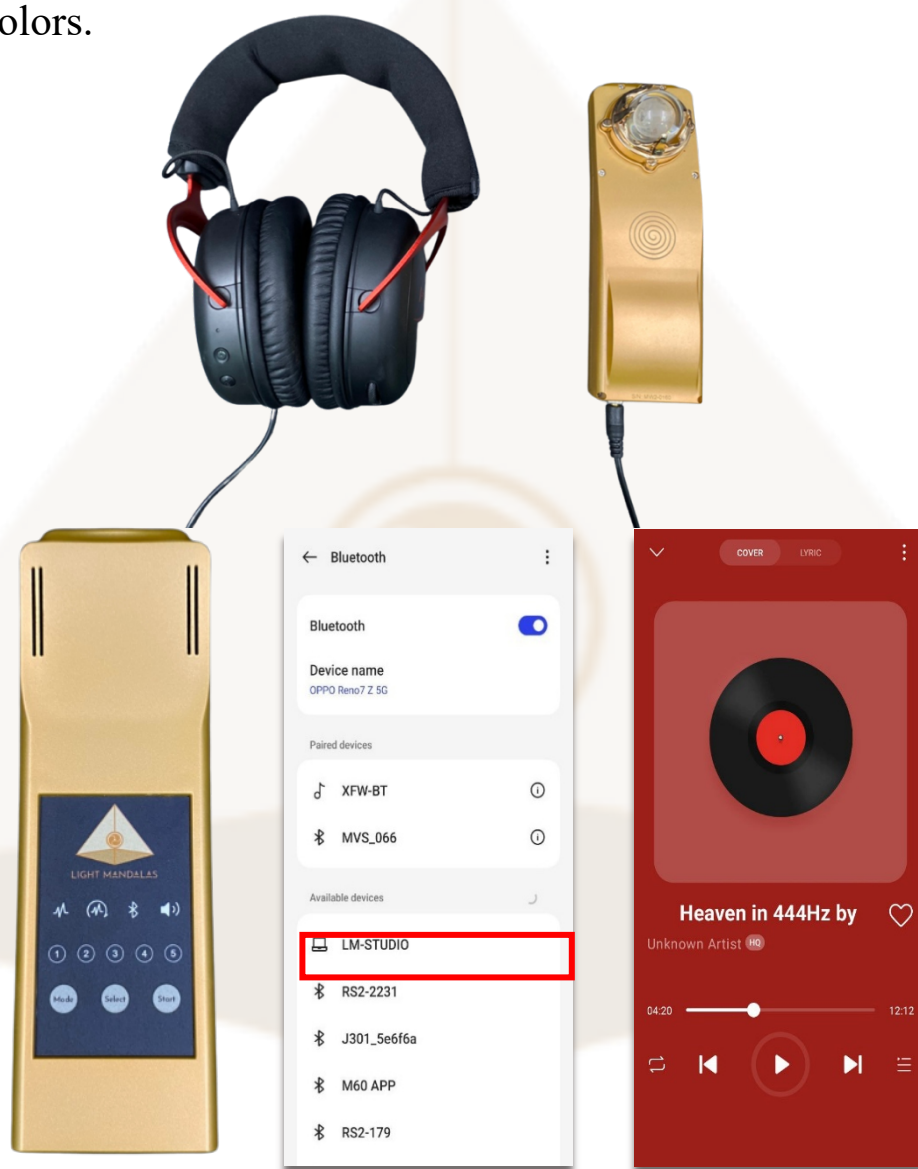

Failures Remedies

Device cannot turn on

- The battery is low. Please charge the device.

Unable to connect with Application

- Make sure the Mandala Wand screen is in Bluetooth Mode.



Disclaimer

Light Mandalas products are designed as energy support tools for relaxation, stress management, and well-being. They are not medical devices and do not replace professional medical treatment or advice.

Any experiences or testimonials shared regarding the use of our products reflect individual responses and are not guaranteed outcomes. Since personal energy and emotional states vary, users may experience different effects, including heightened sensitivity or discomfort. By using our products, you acknowledge that all outcomes are your sole responsibility.

Light Mandalas Co., Ltd. (Thailand) is not liable for any unintended effects, discomfort, or injuries that may arise from product use. If you have any medical concerns, please consult a healthcare professional before using our products.

For any questions, please reach out to us for further support.

Light Mandalas Co., Ltd. (Thailand Headquarters):

- Light Mandalas Website: www.light-mandalas.com
- Email: info@light-mandalas.com
- Instagram: [@mandalawave9](https://www.instagram.com/mandalawave9)
- YouTube Channel: Light Mandalas Official
- TikTok: Light Mandalas
- WeChat: Light Mandalas 光之曼陀罗
- Technical Support Email: cservice@light-mandalas.com